


LOAD CELL TRAINER (STRAIN GAUGE USING LOAD CELL)

Product code: GE 800

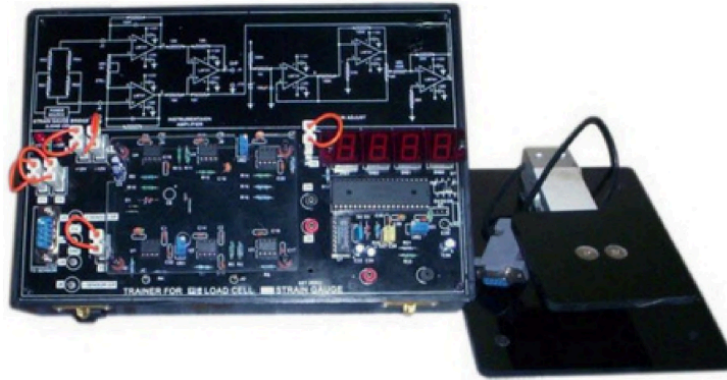

Specifications:

- Study the principle/working of load cell & its applications to measuring Load
- Bridge made up of 4 strain gauges mounted on a specially alloy material is type Load Cell
- Built in excitation voltage for the bridge
- 3.5 digit SSD for direct reading strain in Kg.
- Built in highly sensitive signal amplifier,
- Multi turn potentiometer for zero adjustment,
- Different loads provided for building up of Load.
- Test pins for monitoring at all stages.
- Interconnection: 2/4mm banana sockets
- DC Supply: Built in IC regulated power supplies
- 220V $\pm 10\%$, 50Hz mains operated
- Standard Accessories: Users Manual with patch cords
- Enclosed in an attractive ABS plastic cabinet with cover